# pco.edge 4.2 cooled scMos cameras

lightsheet scanning mode

CLHS FOL USB 3.0

small form factor

high dynamic range 37,500:1

high resolution 2048 x 2048 pixel

> high speed 100 fps

high quantum efficiency up to 82 %

low noise 0.8 electrons









### >> sCMOS image sensor

interfaces >>	CLHS FOL	USB 3.0	
type of sensor	scientific CMOS (sCMOS) monochrome		
resolution (h x v)	2048 x 2048 active pixels		
pixel size (h x v)	6.5 µm x 6.5 µm		
sensor format/diagonal	13.3 mm x 13.3 mm / 18.8 mm		
shutter mode	rolling shutter (RS) with selectable readout modes	rolling shutter (RS) with selectable readout modes global reset - rolling readout (GR)	
MTF	76.9 lp/mm (theoretical)		
fullwell capacity	30,000 e <sup>-</sup>		
readout noise (typ.) <sup>1</sup>	0.8 med e- / 1.3 ms e- @ slow scan 0.9 med e- / 1.4 ms e- @ fast scan	0.8 med e- / 1.3 ms e-	
dynamic range (typ.)	37,500:1 91.5 dB, slow scan	37,500:1 91.5 dB	
quantum efficiency	up to 82 %		
spectral range	300 nm 1100 nm		
dark current (typ.)	< 0.6 e <sup>-</sup> /pixel/s @ 7 °C sensor temperature	< 0.3 e <sup>-</sup> /pixel/s @ 0 °C sensor temperature	
DSNU	< 0.3 ms e <sup>-</sup>	< 0.3 ms e <sup>-</sup>	
PRNU	< 0.3 %	< 0.2 %	
anti blooming factor 2	> 10,000	THE STATE OF THE S	


### >> camera system

interfaces >>	CLHS FOL	USB 3.0	
maximum frame rate @ full resolution	100 fps	40 fps	
exposure/shutter time	100 μs to 10 s (RS)	100 µs to 20 s (RS) 30 µs to 2 s (GR)	
dynamic range A/D <sup>3</sup>	16 bit		
A/D conversion factor	0.46 e-/DN		
pixel scan rate	274.0 MHz fast scan 100.0 MHz slow scan	110.0 MHz	
pixel data rate	548.0 MPixel/s fast scan 200.0 MPixel/s slow scan	220.0 MPixel/s	
binning horizontal	x1, x2, x4		
binning vertical	x1, x2, x4		
region of interest (ROI)	horizontal: steps of 4 pixels vertical: steps of 1 pixel	horizontal: steps of 4 pixels vertical: steps of 1 pixel	
non linearity	< 0.5 %	< 0.6 %	
cooling method	7 °C stabilized, selectable: peltier with forced air (fan) or water cooling (both up to 27 °C ambient)	0 °C stabilized, peltier with forced air (fan) or water cooling (both up to 27 °C ambient)	
trigger input signals	frame trigger, sequence trigger, programmable input (SMA connectors)		
trigger output signals	exposure, busy, line, programmable output (SMA connectors)		
time stamp	in image (1 µs resolution)		

### >>> general

interfaces »	CLHS FOL	USB 3.0	
power delivery	24 VDC (+/- 10 %)		
power consumption	32 W max. (typ. 19 W @ 20 °C)	21 W max. (typ. 12 W @ 20 °C)	
weight <sup>4</sup>	850 g air-cooled	800 g	
	1060 g water-cooled		
operating temperature	+ 10 °C to + 40 °C		
operating humidity range	10 % to 80 % (non-condensing)		
storage temperature range	- 10 °C to + 60 °C		
optical interface	C-mount & F-mount		
lens remote controller	electronic control for Canon EF lenses	not available	
	only air-cooled camera		
maximum cable length	10 km	5 m	
CE/FCC certified	yes		
A CONTRACTOR OF THE PARTY OF TH			

### >> quantum efficiency

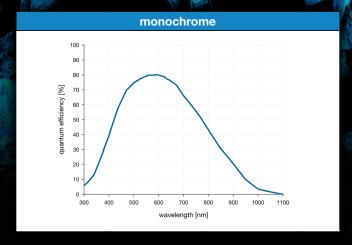

### >> frame rate table<sup>5</sup>

interfaces »	CLHS FOL		USB 3.0
typical examples	fast scan	slow scan	
2048 x 2048	100 fps	35 fps	40 fps
2048 x 1024	200 fps	70 fps	80 fps
2048 x 512	400 fps	140 fps	160 fps
2048 x 256	800 fps	281 fps	315 fps
2048 x 128	1600 fps	562 fps	610 fps
1920 x 1080	189 fps	66 fps	76 fps
1600 x 1200	170 fps	60 fps	69 fps
1280 x 1024	200 fps	70 fps	80 fps
640 x 480	420 fps	150 fps	170 fps
320 x 240	853 fps	300 fps	335 fps


### technical specifications


### pco.edge 4.2

#### dimensions






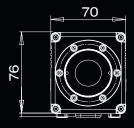


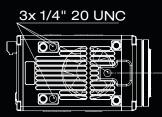

pco.edge USB 3.0



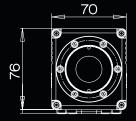

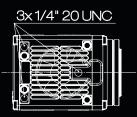

F-mount and C-mount lens adapter are changeable. All dimensions are given in millimeter.



**CLHS FOL** 

air-cooled









The optional Canon lens control adapter enables the user to connect electronic EF- and EF-S Canon lenses allowing to remote control focus and aperture of those lenses.

<sup>1</sup> The readout noise values are given as median (med) and root mean square (rms) values, due to the different noise models, which can be used for evaluation. All values are raw data without any filtering.

2 Based on image sensor data sheet.

The high dynamic signal is simultaneously converted at high and low gain by two 11 bit A/D converters and the two 11 bit values are sophistically merged into one 16 bit value.

4 Measured with C-mount interface.

<sup>&</sup>lt;sup>5</sup> Max. fps with centered ROI.

## technical specifications

### applications



### pco.edge 4.2

bright-field microscopy | fluorescence microscopy | digital pathology | single molecule localization microscopy (SMLM) – PALM, STORM, dSTORM, GSDIM | lightsheet fluorescence microscopy (LSFM) | structured illumination microscopy (SIM) | calcium imaging | förster resonance energy transfer (FRET) | fluorescence recovery after photobleaching (FRAP) | high-speed bright-field ratio imaging | high throughput screening | high content screening | biochip reading | total internal reflection microscopy (TITF) | spinning disk confocal microscopy | 3D metrology | ophthalmology | photovoltaic inspection | industrial quality inspection | wafer inspection | image intensifier imaging | lucky astronomy | particle tracking velocimetry (PTV)



With pco.camware you control all camera settings, the image acquisition and the storage of your image data. The pco.sdk is the complementary software development kit. It includes dynamic link libraries for user customization and integration on Windows-PC platforms. Drivers for popular third party software packages are also available for you.

All this items like pco.camware, pco.sdk and third party drivers, are free-to-download at

### third party integrations



### contact

**pco europe** +49 9441 2005 50 info@pco.de pco.de

#### pco america

+1 866 678 4566 info@pco-tech.com pco-tech.com

#### pco asia

+65 6549 7054 info@pco-imaging.com pco-imaging.com

### pco china

+86 512 67634643 info@pco.cn pco.cn



for application stories please visit our website



An Excelitas Technologies Brand